# Angelcare®

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**AC601**

**OWNER'S MANUAL**  
Keep for future reference

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KEEP THIS INSTRUCTION  
MANUAL FOR FUTURE REFERENCE

FAILURE TO FOLLOW THE WARNINGS  
AND THE INSTRUCTIONS COULD RESULT  
IN SERIOUS INJURY OR DEATH  
OR DAMAGES TO THE PRODUCT

# ⚠ WARNING

## IMPORTANT SAFETY INSTRUCTIONS

THIS PRODUCT CANNOT REPLACE RESPONSIBLE ADULT SUPERVISION. This product is designed to aid in the monitoring of your child. When using the monitor, you must still provide the supervision necessary for the continued safety of your child. Monitoring of premature babies, or those considered to be at risk, should always be carried out under the supervision of a health professional. This monitor is not intended to be used as a medical monitor.

**STRANGULATION HAZARD.** Children have STRANGLED in cords. Keep the adapter cords and sensor pad cord out of the reach of children. Do not put the Nursery Unit inside a crib or play yard. Because of the serious strangulation risk to infants and toddlers, parents and caregivers should never place the nursery unit within 3 feet (1 meter) of a crib. NEVER use extension cords with AC adapters. ONLY use the AC adapters provided.

Route sensor pad cord only below crib mattress and out of reach of children.

When an AC adapter is plugged into an electrical outlet, do NOT touch the exposed end.

ALWAYS be sure that both the transmitter and receiver are working properly and are within range of each other. Test monitor before first use, periodically and when changing location of the camera.

ALWAYS be sure to immediately check on baby whenever an alarm sounds.

Do NOT use monitor near water (such as bath tub, sink, ect.).

Keep monitor AWAY from heat sources (such as stoves, radiators, ect.).

ALWAYS place both the Nursery Unit and Parent Unit on a flat surface in an upright position and allow for proper ventilation. Do NOT place on sofas, cushions, beds, etc. which may block ventilation.

This monitor uses public airwaves to transmit signals. The Parent Unit may pick up signals or interference from other nursery monitors in your area. Your monitor may even be picked up by other households. To protect your privacy, make sure both units are turned off when not in use. If you experience interference, change the channel on both units.

This product is NOT a toy. Do NOT allow children to play with it.

This product contains small parts. Adult assembly required. Exercise care when unpacking and assembling the product and keep small parts away from children.

Use and store the monitor at a temperature between 50°F (10°C) and 104°F (40°C).

## BATTERY WARNING:

Use ONLY AAA (LR03) RECHARGEABLE batteries for the Parent Units.

Do NOT mix old and new batteries.

Do NOT mix alkaline, standard or rechargeable batteries.

Insert batteries with correct polarity.

Do NOT short-circuit supply terminals.

ONLY use the recommended batteries or equivalent (volts and size) and the correct adapter.

Remove batteries when product is stored for long periods of time or when batteries are exhausted.

Dispose exhausted batteries properly.

Do NOT recharge non-rechargeable batteries.

Rechargeable batteries are ONLY to be charged under adult supervision.

Keep all batteries away from children.

## IMPORTANT SAFETY INSTRUCTIONS

1) Read these instructions- All the safety instructions should be read before the product is operated.

2) Keep these instructions. The safety and operating instructions should be retained for future reference.

3) Heed all warnings – All warnings on the product and in the operating instructions should be adhered to.

4) Follow all instructions – All operating and use instructions should be followed.

5) Water and moisture- Do NOT use this apparatus near water. For example, near a bath tub, wash bowl, kitchen sink, laundry tub, in a wet environment or near a swimming pool and the like. Care should be taken so that objects do not fall onto units and that liquids are not spilled into either unit or sensor pad through its openings.

6) Cleaning - Disconnect all units before cleaning. Do NOT immerse in water any part of the Nursery Unit, Parent Units, or docking station. Clean only with dry cloth. Do NOT spray with cleaning products or solvents. Wipe sensor pad with a cloth moistened very slightly with an antiseptic or mild detergent.

7) Ventilation - Do NOT block any ventilation openings. Allow for proper ventilation when units are in use. Do NOT cover the Nursery Unit or Parent Units with any object such as a blanket. Do NOT place it in a drawer or in any location which would muffle the sound or interfere with the normal flow of air.

8) Heat sources- Do NOT install near any heat sources such as radiators, heat registers, stoves, TVs or other apparatus (including amplifiers) that produce heat. Heat can damage the case or electrical parts.

9) Power cord protection- Protect the sensor pad cord and adapter's power cord from being walked on or pinched by furniture or other items – Power supply cords should be routed so that they are not likely to be walked on or pinched by items upon or against them, paying particular attention to cords and plugs, convenience receptacle, and the point where they exit from the product.

10) Accessories- Only use attachments/accessories specified by the manufacturer. Do NOT place this product on an unstable cart, stand, tripod, bracket or table. The product may fall, causing serious injury to a child or adult and serious damage to the product. Use only with a cart, stand, tripod, bracket or table recommended by the manufacturer or sold with the product. Any mounting of the product should follow

the manufacturer's instructions.

11) Unplug this appliance during lightning storms or when not used for long periods of time.

12) Servicing- Do NOT attempt to service this product yourself as opening or removing covers (except battery covers) may expose you to dangerous voltage or other hazards.

13) Damage required servicing- Unplug this product from the wall outlet and refer all servicing to qualified service personnel. Servicing is required when the apparatus has been damaged in anyway, such as power supply cord or plug is damaged, liquid has been spilled or objects have fallen into the apparatus, the apparatus has been exposed to rain or moisture, does not operate normally, or has been dropped.

14) Replacement parts - When replacement parts are required, be sure the service technician has used replacement parts specified by the manufacturer or have the same characteristic as the original part. Unauthorized substitution may result in fire, electric shock or other hazards.

15) Safety check- Upon completion of any service or repair to the product, ask the service technician to perform safety checks to determine that the product is in proper operating condition.

16) Usage-This baby monitor is intended as an aid. It is NOT a substitute for responsible and proper adult supervision and should not be used as such. Use this appliance ONLY for its intended use as described in this manual.

17) Power Sources- Use ONLY with 110-120 volt AC electrical outlets.

18) Overloading- Do NOT overload wall outlets, extension cords or integral convenience receptacles as this can result in a risk of fire or electric shock.

19) Object and Liquid Entry- NEVER push objects of any kind into this product through casing slots as they may touch dangerous voltage points or short out parts that could result in a risk of fire or electric shock. NEVER spill liquid of any kind on the product.

20) Antenna -The antenna used for this transmitter must be installed to provide a separation distance of at least 20 cm (7.9") from all persons and must not be collocated or operating in conjunction with any other antenna or transmitter.

21) Power Lines- An outside antenna system should NOT be located in the vicinity of overhead power lines or other electric light or power circuit, or where it can fall into such power lines or circuits as contact with them may be fatal.

22) STRANGULATION HAZARD. Children have STRANGLED in cords. Keep the adapter cords and sensor pad cord out of the reach of children. Because of the serious strangulation risk to infants and toddlers, parents and caregivers should NEVER place the nursery unit within 3 feet (1 meter) of a crib. NEVER use extension cords with AC adapters. ONLY use the AC adapters provided.

23) Interference- This monitor uses public airwaves to transmit signals. The Parent Unit may pick up signals or interference from other nursery monitors in your area. Your monitor may even be picked up by other households. To protect your privacy, make sure both units are turned off when not in use. If you experience interference, change the channel on both units.

24) Vibration Source- When the monitor is in use, do NOT use a crib mobile or other crib accessories that may produce vibrations. Other sources of vibrations may include fans, washing machines, loud music, etc. These sources of vibration may interfere with the performance of your Angelcare® monitor and the alarm will not sound if it continues to detect any movement other than your baby's. Be sure all vibration sources are eliminated before using the monitor.

You may experience false alarms for several reasons. The most common reason is the Nursery Unit has not been turned off after you have removed your baby from the crib. Another reason is the Sensor Pad has not been placed under the mattress properly.

Be sure the Sensor Pad is on a flat, firm surface with

the printed side facing up. The sensitivity level may also need to be adjusted. Refer to the Movement Sensitivity Adjustment section on P 9. Always be sure to immediately check on baby whenever an alarm sounds.

25) Sensor Pad- Your sensor pad can be used in any crib that is stationary and designed to safely hold a baby. If the crib has a spring or slat base, place a hard board between the crib base and the Sensor Pad measuring at least 13" x 13" (33 x 33 cm) and 1/4 inch (6 mm) thick. Only use your Angelcare® monitor in a stationary sleeping environment such as a crib where the Sensor Pad can rest on a completely flat, firm surface.

ALWAYS route the sensor pad cord below crib mattress ONLY and out of reach of children.

## IMPORTANT SAFETY INSTRUCTIONS RELATED TO USE OF BATTERIES

1) Use ONLY the type and size of batteries specified in the user manual. There is a risk of explosion if batteries are replaced with an incorrect type.

2) Use ONLY AAA (LR03) rechargeable batteries and/or the AC adapters provided for the Parent Units. You must use only rechargeable batteries in the Parent Unit. DO NOT use alkaline batteries or you will severely damage the unit.

3) Charge the batteries for 15 hours prior to first use. Be sure to have the monitor in the OFF position during this period. ONLY use the AC adapters included with your monitor.

4) Do NOT use these adapters with any other equipments.

5) When unplugging the AC adapter from the nursery unit, be sure it is switched to the off position. If it isn't, the unit will revert to battery mode if batteries are installed.

6) Do NOT use rechargeable batteries in the Nursery Unit.

7) Do NOT dispose batteries in fire. They may explode. Check with local regulations for possible special disposal instructions.

8) Exercise care in handling batteries in order not to short out the batteries with conducting materials such as rings, bracelets and keys. Overcharging, short circuiting, reverse charging, mutilation, or incineration of batteries must be avoided to prevent one or more of the following occurrences: release of toxic materials, release of hydrogen and/or oxygen gas, rise in surface temperature.

9) Do NOT attempt to recharge batteries provided with or identified for use with this product that are not intended to be charged. The batteries may leak corrosive electrolyte or explode.

10) Do NOT attempt to rejuvenate the batteries provided with or identified for use with this product by heating them. Sudden release of battery electrolyte may occur, causing burns or irritation to eyes or skin.

11) If and when necessary, the batteries must be replaced with identical new ones or equivalent (volts and size). If a battery to be replaced is connected with other batteries in series, replace the other batteries at the same time as well. Do NOT mix old and new batteries.

12) Reverse polarity installation of the batteries in the end product must be avoided. Reverse insertion of batteries can cause charging, and that may result in leakage or explosion.

13) Remove batteries from this product if you store the product over 30 days because the batteries could leak and damage the product.

14) Discard 'dead' batteries as soon as possible since 'dead' batteries are more likely to leak in a product. If a battery has leaked or vented, it should be replaced immediately using protective gloves.

15) Do NOT store this product, or the batteries provided with or identified for use with this product, in high-temperature areas. Batteries that are stored in a freezer or refrigerator for the purpose of extending shelf life should be protected from condensation during storage and defrosting. Batteries should be stabilized at room temperature prior to use after cold storage.

## WARNING

Modifications not authorized by the manufacturer may void users authority to operate this device.

**NOTE:** This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:


- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) This device must accept any interference received, including interference that may cause undesired operation.

## Features


### Parent Unit
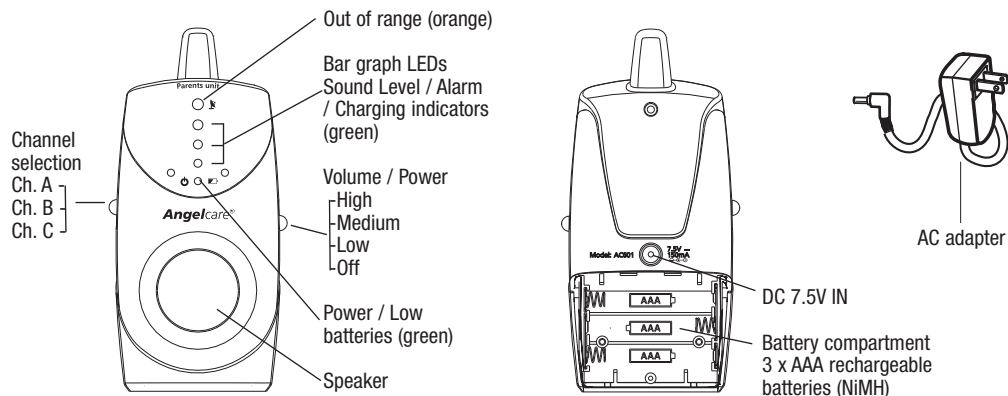
### Installing Batteries in Parent Unit

- The Parent Unit requires three AAA rechargeable NiMH batteries (included). **Note:** Use only AAA rechargeable NiMH batteries in the Parent Unit. DO NOT use regular alkaline batteries, they may explode.
- To insert batteries, remove the battery cover on the back of the unit and insert the batteries matching positive (+) to positive (+) and negative (-) to negative (-) and replace the battery cover.
- Leave the AC adapter plugged to the Parent Unit for 15 continuous hours before using your monitor to fully charge rechargeable batteries.
- Recharge the batteries when the Low Battery Indicator flashes by plugging in the AC adapter.
- The AC adapter can also be used as an alternative operation power source.
- IMPORTANT:** Use only the AC adapters included with your monitor. Do not use these adapters with any other equipment.

## Nursery Unit
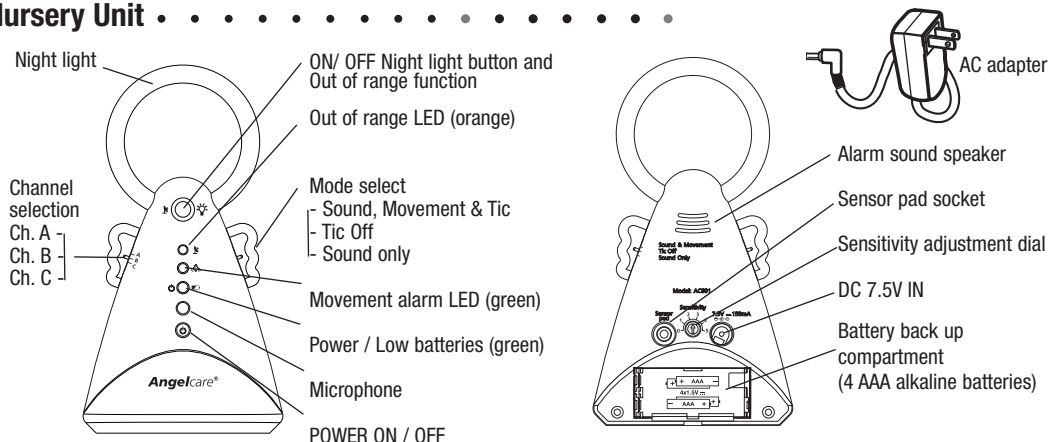
### Installing Batteries in Nursery Unit

- The Nursery Unit requires four AAA alkaline batteries in case of AC power failure (included).
- To insert batteries, remove the battery cover on the bottom of the unit and insert the batteries matching positive (+) to positive (+) and negative (-) to negative (-) and replace the battery cover.
- Insert the AC adapter into the Nursery Unit and plug the adapter into an electrical outlet. Unit will run off of the AC adapter unless there is a power outage, in which case the unit will begin using the batteries.
- Replace the batteries when the Low Battery Indicator flashes. If the batteries are not replaced, the monitor will stop working.
- IMPORTANT:** Use only the AC adapters included with your monitor. Do not use these adapters with any other equipment.

# Setting Up Your Monitor . . . . .

## SENSOR PAD (See figs. 1 and 2)

Place the Sensor Pad under the center of the crib mattress with the printed side facing up. **Make sure there is NO bedding between the mattress and the Sensor Pad. The Sensor Pad must rest on a completely flat and firm surface. If the crib has a spring or slat base, place a hard board between the crib base and the Sensor Pad measuring at least 13" x 13" (33x33 cm) and 1/4 inch (6 mm) thick.** If your crib has a plywood base, your mattress support is firm enough to use the monitor properly. Only use your Angelcare® monitor in a stationary sleeping environment such as a crib where the Sensor Pad can rest on a flat, firm surface. Be sure to place the printed side of the Sensor Pad facing up. Now plug the Sensor Pad cord into the jack on the back of the Nursery Unit marked "Sensor Pad". Secure the sensor pad cord to crib using a tie wrap. Run the cord on the floor to ensure it is out of the baby's reach. Turn the Nursery Unit on using the switch on the left. The green Power Indicator will light up. NOTE : The Sensor Pad works with all baby bed sizes and all mattresses, except memory foam mattresses or mattresses with a hollow frame. You will need to test your Sensor Pad after the set-up is completed. Once you are done using the Sensor Pad, please remove it from the baby's crib and store it away. The Angelcare movement monitor is not designed for use with adults, children or infants older than 24 months.


## ! WARNING

**STRANGULATION HAZARD - Children have STRANGLED in cords. Keep the adapter cords and sensor pad cord out of the reach of children. Do not put the Nursery Unit inside a crib or play yard. Because of the serious strangulation risk to infants and toddlers, parents and caregivers should never place the nursery unit within 3 feet (1 meter) of a crib. NEVER use extension cords with AC adapters. ONLY use the AC adapters provided.**

Figure 1



## PLACEMENT OF SENSOR PAD CORD

### CORRECT

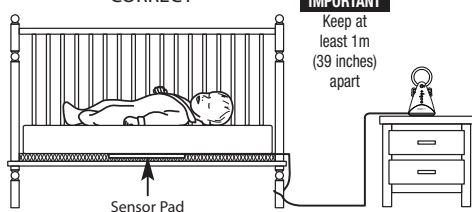**IMPORTANT**  
Keep at  
least 1m  
(39 inches)  
apart

### INCORRECT

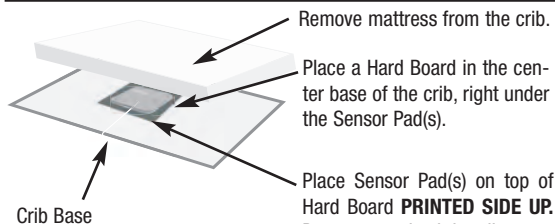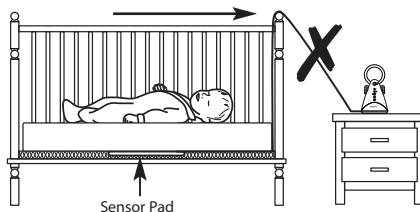
Remove mattress from the crib.

Place a Hard Board in the center base of the crib, right under the Sensor Pad(s).


Place Sensor Pad(s) on top of Hard Board **PRINTED SIDE UP**. Put mattress back in crib.

### 1 Sensor Pad

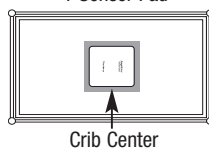Crib Center

### 2 Sensor Pads

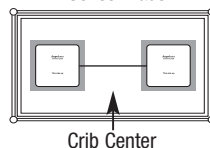Crib Center

**The Sensor Pad works with most cribs and mattresses, except memory foam and hollow frame mattresses.**

## SELECT YOUR DETECTION MODE

On the right of the Nursery Unit select one of the following three functions by moving the wing upwards or downwards:

- **Sound & Movement:** You will be able to hear your baby's sounds and see sound lights. Your baby's movements are monitored and you will hear a "Tic" with every movement your baby makes. If no movement is detected, the "Tic" sound will stop and after 20 seconds the alarm will sound.
- **Tic Off:** In this mode you will still hear your baby's sounds and see sound lights, and your baby's movement is still being detected, but no "Tic" sound will be heard. If no movement is detected after 20 seconds the alarm will sound.
- **Sound Only:** The Movement Sensor is disabled but you will still be able to hear your baby's sounds and see sound lights. Your Angelcare® works like a regular sound monitor when set to "Sound Only".

## Testing Your Monitor . . . . .

**Always test your monitor before using it the first time and then on a periodic basis. Before testing your monitor, complete the setup process as described above.**

### TESTING THE MOVEMENT FUNCTION

Note: Make Sure your Sensor Pad is set up before proceeding with this section. If you haven't set up your Sensor Pad please go to page 5 for instructions on how to properly set up the Sensor Pad.

Make sure the Sensor Pad's cord is plugged into the Nursery Unit. ***Switch the Nursery Unit to the "Sound / Movement" function.*** Turn the Parent Unit ON and select either channel A, B or C. Make sure it is set to the same channel as the Nursery Unit. The green power indicator light will light up. Gently pat your hand on the mattress. The green Movement Indicator light on the Nursery Unit will flash continuously.

Remove your hand from the mattress. Due to the lack of movement, the Movement Indicator light will stop flashing, and a pre-alarm beep will sound after 15 seconds. Five seconds later the alarm will sound, the green Alarm Indicator Light on the Nursery Unit will flash continuously and the green lights on the Parent Unit will also flash continuously.

**IMPORTANT:** If the alarm does not sound when you are not touching the mattress, the Sensor Pad may be detecting a continuous vibration from another source from inside and/or outside your baby's room. Common sources of this interference are crib mobiles, fans, washing machines or even loud music. Therefore, it is very important to make sure all vibration sources are eliminated before using the monitor.

Place the crib near a solid supporting wall and away from strong air currents. Do not use crib mobiles or accessories that may cause vibrations. Avoid contact with the crib when the monitor is in use. If you find the alarm will still not sound when you are not touching the mattress or crib, the monitor is on, and you have removed all external sources of vibration, the sensitivity level may need to be decreased (Refer to the movement sensitivity adjustment section on page 7).

Test the stopping of the alarm: to stop the alarm, place your hand lightly on the mattress so that the Sensor Pad starts detecting movement again. You may also shut the Nursery Unit OFF to stop the alarm.

### TESTING THE "TIC" FUNCTION

Make sure the Sensor Pad's cord is plugged into the Nursery Unit. ***Switch the Nursery Unit to the "Sound / Movement" function.*** Turn the Parent Unit ON and select either channel A, B or C. Make sure it is set to the same channel as the Nursery Unit. The green power indicator light will light up. Gently pat your hand on the mattress. The green Movement Indicator light on the Nursery Unit will flash and you will hear the Tic sound with every movement.

Remove your hand from the mattress. Due to the lack of movement, the ticking will stop, the Movement Indicator light will stop flashing and a pre-alarm beep will sound after 15 seconds. Five seconds later the alarm will sound, the green alarm indicator light on the Nursery Unit will flash continuously and the green lights on the Parent Unit will also flash continuously. To test the sound function, please refer to the Sound Only section below.

***The "Tic" sound will be heard when the Sound & Movement option is chosen. If you do not wish to hear the "tic", select "tic off" on the right wing of the Nursery Unit.***

**TESTING THE SOUND FUNCTION**


**Switch the Nursery Unit to "Sound Only".** Turn the Parent Unit ON and select either channel A, B or C. Make sure it is set to the same channel as the Nursery Unit. The green power indicator LED will light up. Choose the channel that gives you the best reception.

Test the sound: have someone else speak into the Nursery Unit or place a radio near the Nursery Unit. You should hear sound from the Parent Unit. If there is no sound, refer to the Troubleshooting section. Adjust the volume to the desired level using the slide button on the right of the Parent Unit.

**TESTING THE BATTERY POWER**

Verify that the units are operating on battery power by installing 3 AAA NiHM rechargeable batteries (included) in the Parent Unit and 4 regular AAA batteries (included) in the Nursery Unit and unplugging the AC adapters. The green Power Indicator should light up on all units.

**Movement Sensitivity adjustment • • • • •**

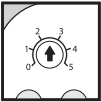**Movement Sensitivity Adjustment**

The movement sensitivity adjustment dial has been preset to a level which will suit most situations. The highest sensitivity level is 5 (most sensitive) and the lowest sensitivity level is 0 (least sensitive).

You may get false alarms if your baby is in a very deep sleep or has moved to the edge of the crib. If you get false alarms due to these circumstances, you should increase the sensitivity level by using the sensitivity dial on the back of the Nursery Unit.

**NOTE:** Adjust the sensitivity level if you receive false alarms and you have followed the advice in the Troubleshooting section for false alarms. Make a note of the position of the sensitivity dial before you change it should you wish to return to that setting.

To increase the sensitivity level, start by increasing it by half a level. For example, if the level is currently at 2, increase it to 3. Use the monitor as you normally would. If you still experience false alarms, increase it by another level. If you are still experiencing false alarms when you have reached level 5, please consult our customer service link on our website at [www.angelcare-monitor.com](http://www.angelcare-monitor.com)

Because the Sensor Pad may detect movement from many sources from both inside and outside your baby's room, be sure all vibration sources are eliminated before changing the movement sensitivity level. Check to see if the Sensor Pad is resting on a flat, firm surface. Place the crib near a solid supporting wall and away from strong air currents. Do not use crib mobiles or accessories that may cause vibrations. Avoid contact with the crib when monitor is in use.

To make sure you have avoided all external sources of vibration, take your baby out of the crib while the monitor is on. If the alarm does not sound after you take your baby out of the crib, the sensitivity level may need to be decreased. To decrease the sensitivity level, adjust by a level by going from 2 to 1. If you are still experiencing this problem when the level is set to 0, please consult our customer service link on our website at [www.angelcare-monitor.com](http://www.angelcare-monitor.com)

## Special features • • • • • • • • • •

### Out of range indicator

Your monitor is equipped with an optional “Out of range” indicator function, which will advise you when the two units are not communicating with each other. Out of range situations can be caused by distance or by physical interferences.

When the units are out of range, you will hear a DOUBLE BEEP every 8 seconds and the “out of range” icon will flash on the Parent Unit.

Range: up to 820 ft (250 meters).

To activate this function, make sure that both units are ON. Press and hold the night light button on the Nursery Unit for 2 seconds. The orange lights next to the out of range icons on both units will turn orange to remind you that the function has been activated.

### IMPORTANT:

- If you turn the Nursery Unit off, you will have to reactivate this option when you turn it on again.
- If you turn the Parent Unit off, you will have to reset and reactivate this option when you turn it on again. To reset it, press and hold the night light button on the Nursery Unit for 2 seconds. Then press and hold the night light button on the Nursery Unit again for 2 seconds to activate the function.

**Note:** If you turn the Nursery Unit off when the out of range option is activated, the out of range alarm will beep as the Parent Unit cannot detect the Nursery Unit anymore. To stop it, you will have to either turn the Nursery Unit ON or turn the Parent Unit off.

### Sound Lights

The Parent Unit features Sound Lights so you can “see” your baby’s sounds if background noise prevents you from hearing your baby. As the sounds get louder, more lights illuminate.

### Movement Indicator light

The Movement Indicator Light is located only on the Nursery Unit. When the Nursery Unit is set to “Sound & Movement” or to “Tic Off”, the Movement Indicator Light will flash green with your baby’s every movement. The flashing green light will let you know that movement is being detected. This will allow you to check on your sleeping baby in the dark without disturbing.

### Night light

The halo on the Nursery Unit acts as a soothing nightlight and allows you to check your baby in the dark with minimal disruption. Press the halo button to turn the night light ON.

## Tips For Using Your Monitor • • • • •

- When the monitor is in use, do not use a crib mobile or other crib accessories that may produce vibrations. Other sources of vibrations may include fans, washing machines, loud music, etc. These sources of vibration may interfere with the performance of your Angelcare® monitor and the alarm will not sound if it continues to detect any movement. Be sure all vibration sources are eliminated before using the monitor.
- There are several reasons you may experience false alarms. The most common reason is the Nursery Unit has not been turned off after you have removed your baby from the crib. Another reason is the Sensor Pad has not been placed under the mattress properly. Be sure the Sensor Pad is on a flat, hard surface with the printed side facing up. The sensitivity level may also need to be adjusted. Refer to the Movement Sensitivity Adjustment section.
- Your monitor can be used in any crib that is stationary and designed to safely hold a baby. Make sure there is no bedding between the mattress and the Sensor Pad. The Sensor Pad must rest on a completely flat and firm surface. If the crib has a spring or slat base, place a hard board between the crib base and the Sensor Pad measuring at least 13" x 13" (33x33 cm) and 1/4 inch (6 mm) thick. Only use your Angelcare® monitor in a stationary sleeping environment such as a crib where the Sensor Pad can rest on a hard, flat surface.
- Please note the difference between the Low Battery Indicator and the Alarm Indicator. The low battery indicator is a rapid ticking sound on the Parent Unit and the green low battery indicator light slowly flashes. The Alarm Indicator is a loud beeping sound and the green alarm indicator light flashes rapidly.
- When unplugging the AC adapter from each unit, be sure the unit is switched to the off position. If it isn't, the unit will revert to battery mode if batteries are installed.
- You can use rechargeable batteries in the Nursery Unit, however, a rechargeable battery cannot be recharged by plugging in the AC adapter. Remove batteries from the unit to recharge.
- The range of this monitor is approximately 820 feet (250 meters). However, the range will vary depending on the surroundings. Test the range before use.

## Troubleshooting • • • • •

Problem	Possible Cause	Solution
False alarms	<ul style="list-style-type: none"> <li>Baby was removed from the crib and the Nursery Unit was not switched off.</li> <li>The Sensor Pad is not detecting your baby's movements because the sensitivity level is too low.</li> <li>Sensor Pad is not in full contact with the crib mattress or is not resting on a flat, firm surface.</li> <li>Baby has moved far away from the Sensor Pad.</li> <li>Sensor Pad cord is not connected properly to the Nursery Unit.</li> </ul>	<ul style="list-style-type: none"> <li>Turn the Nursery Unit off.</li> <li>The sensitivity level may need to be increased. Refer to the Movement Sensitivity Adjustment section on Page 7.</li> <li>Make sure there is no bedding between the mattress and the Sensor Pad. The Sensor Pad must rest on a completely flat and firm surface. If the crib has a spring or slat base, place a hard board between the crib base and the Sensor Pad measuring at least 13" x 13" (33x33 cm) and 1/4 inch (6 mm) thick."</li> <li>Refer to the Movement Sensitivity Adjustment section on Page 7.</li> <li>Unplug and reconnect the cord.</li> </ul>
Alarm will not sound	<ul style="list-style-type: none"> <li>Nursery Unit is detecting movement when crib is being touched.</li> <li>Nursery Unit is detecting movement from outside the crib such as motorized devices, strong drafts, etc.</li> <li>Sensor Pad sensitivity level is set too high.</li> </ul>	<ul style="list-style-type: none"> <li>Avoid contact with the crib when monitor is on.</li> <li>Place crib near a solid supporting wall to avoid Sensor Pad from detecting movement from outside the crib. Turn off crib mobile, fans, loud music or anything else that may cause vibrations. You may need to decrease the sensitivity.</li> <li>Refer to the Movement Sensitivity Adjustment section on Page 7.</li> </ul>
No sound is being transmitted.	<ul style="list-style-type: none"> <li>Parent and Nursery Units are too far apart.</li> <li>Batteries are low or improperly installed.</li> <li>AC adapters are not connected properly or electrical outlet is out of order.</li> <li>One of the units is not on.</li> <li>Units are set to different channels.</li> <li>Volume level is too low</li> </ul>	<ul style="list-style-type: none"> <li>Bring the Parent Unit closer to Nursery Unit.</li> <li>Check or replace batteries.</li> <li>Check connections or change to another electrical outlet.</li> <li>Make sure both units are on.</li> <li>Make sure both units are on the same channel, either on A, B or C.</li> <li>Increase volume level.</li> </ul>
Power Indicator (green light) does not light up.	<ul style="list-style-type: none"> <li>Batteries are low or improperly installed.</li> <li>AC adapters are not connected properly or electrical outlet is out of order.</li> <li>One of the units is not on.</li> </ul>	<ul style="list-style-type: none"> <li>Check or replace batteries.</li> <li>Check connections or change to another electrical outlet.</li> <li>Make sure both units are on the same channel, either on A, B or C.</li> </ul>
Very weak reception.	<ul style="list-style-type: none"> <li>Parent Unit and Nursery Unit are too far apart.</li> <li>One or both units are placed near a large metal mass.</li> <li>One or both units are not in the upright position.</li> <li>Batteries are low.</li> </ul>	<ul style="list-style-type: none"> <li>Bring the Parent Unit closer to Nursery Unit.</li> <li>Change the position of one or both units.</li> <li>Place units on flat surface away from baby's reach.</li> <li>Check or replace batteries.</li> </ul>

Problem	Possible Cause	Solution
There is static, distortion, interference, etc	<ul style="list-style-type: none"><li>• Parent Unit and Nursery Unit are too far apart.</li><li>• Units are not set to the same channel.</li><li>• One of the units is not on.</li><li>• The Parent Unit is located near motorized devices, fluorescent lights, television, etc.</li><li>• Batteries are low.</li></ul>	<ul style="list-style-type: none"><li>• Bring the Parent Unit closer to Nursery Unit.</li><li>• Make sure both units are on the same channel, either on A, B or C.</li><li>• Make sure both units are on.</li><li>• Rotate Parent Unit or move it away from source of interference.</li><li>• Recharge the batteries.</li></ul>
There is feedback (squealing noises).	<ul style="list-style-type: none"><li>• Units are too close to each other.</li><li>• The volume is too high.</li></ul>	<ul style="list-style-type: none"><li>• Move units farther apart (at least 10 feet apart).</li><li>• Decrease volume on Parent Unit.</li></ul>

Limited Warranty . . . . .

Angelcare® warrants to the original owner of this product that this Angelcare®monitor (model AC601) is free from any material and workmanship defects, for a period of one (1) year from the date of the original purchase. If this monitoring system fails to function properly when used as directed under normal conditions within the one-year period, Angelcare® will either repair or replace the product, at our discretion, free of charge. Please contact your distributor.

The product must be accompanied by a proof of purchase, either a bill of sale or other proof that the monitoring system is within the warranty period. Angelcare® will bear the cost of repairing or replacing the product and shipping it back to you.

**This warranty does not apply to a product which has been damaged as a result of improper maintenance, an accident, improper voltage supply or any other form of misuse. The warranty is also void if the owner repairs or modifies the product in any way. Angelcare® is not liable for any incidental or consequential damages with regard to this product. The warranty also excludes any liability other than what is stated above. No other warranty is given.**

LEGAL RIGHTS VARY FROM COUNTRY TO COUNTRY. SOME COUNTRIES DO NOT ALLOW THE EXCLUSION OR LIMITATION OF INCIDENTAL OR CONSEQUENTIAL DAMAGES, SO THE ABOVE LIMITATIONS MAY NOT APPLY TO YOU.

WE ARE NOT IN ANY WAY RESPONSIBLE FOR ANY OUTCOME RESULTING FROM THE USE OF THIS PRODUCT.

REGISTER ANGELCARE PRODUCTS

Register your Angelcare monitor and other Angelcare products on our website: [www.angelcare-monitor.com](http://www.angelcare-monitor.com)

3980, rue Saint Ambroise  
Montréal, Québec, H4C 2C7,  
Canada

**IMPORTANT! This Angelcare® monitor is a personal care product. Please do not return this product to the store. If you have questions or need help, please call Angelcare before returning this product.**

Replacement Parts and Warranty Information (USA & Canada only)  
[www.angelcare-monitor.com](http://www.angelcare-monitor.com) - [angelcare@cdecinc.ca](mailto:angelcare@cdecinc.ca)  
**1-877-281-8827**

**Protect the environment!**  
This electronic equipment contains reusable material and must not be disposed with domestic waste. Please take it to a collecting point for electric and electronic equipment (ask your local authorities for details).

**Power supply:** Four AAA alkaline batteries (not included) or 7.5V AC adapter for Nursery unit. Three AAA rechargeable batteries and 7.5V AC adapter for the Parent Unit.

**CAUTION:** To reduce the risk of electric shock, do not expose this product to rain or moisture.

Danger of explosion if batteries are incorrectly replaced. Batteries included: AAA NIMH rechargeable batteries, manufactured by GPI INTERNATIONAL LIMITED (GN60AAAHC, 1.2V 600mAh) and by HIGHPOWER TECHNOLOGY CO., LTD (HFR-44AAAJ600, 1.2V 600mAh). **They may only be replaced by AAA NIMH rechargeable batteries.**

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